


MEISINGER

Easy-Clean Scissors

Surgical Scissors made of detachable individual parts –
for easy and safe processing

The MEISINGER Easy-Clean Scissors made of high-quality stainless steel can be easily and completely disassembled and put back together again. Residues below the hinge can be completely removed – for a residue-free processing*. For dismantling, the scissors only need to be opened at a right angle. Then the upper part of the scissors is lifted briefly. In addition, the scissors have a black coating that is applied using PVD procedure (vacuum-based coating). Thanks to this black coating the products reach high hardness and durability. Furthermore, disturbing light reflections are reduced, and glare-free work is possible.

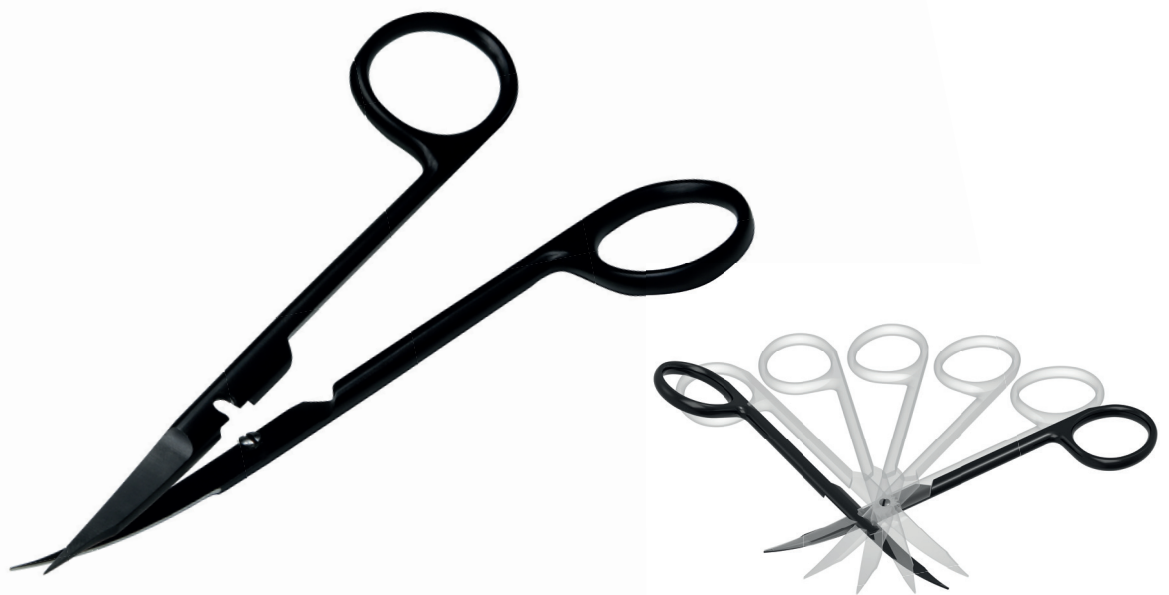

At a glance

- Surgical scissors made of high-quality stainless steel
- Easy and complete dismantling for residue-free processing*
- Black coating for high hardness and durability as well as glare-free work
- Scissors in two different versions





* In correspondence with required processing procedures according to EN ISO 17664.

Our sales force would be pleased to
answer your questions: +49 2131 2012-303

Meisinger
since
1888

MEISINGER Easy-Clean Scissors

Fax for ordering +49 2131 2012-222

	Figure	Length	Pieces / Kit	Amount 
	CM01S	-	-	
 illustrated 1:2	CM001	12 cm	1	
 illustrated 1:2	CM002	13 cm	1	

CE 0044

Available from the dental trade.

Your trade depot:

Your delivery address:

Hager & Meisinger GmbH
 Hansemannstr. 10 • 41468 Neuss • Germany
 Phone: +49 2131 2012-0 • Fax: +49 2131 2012-222
 info@meisinger.de • www.meisinger.de

Meisinger USA, L.L.C.
 10150 E. Easter Avenue • Centennial, Colorado 80112 • USA
 Phone: +1 (303) 268-5400 • Fax: +1 (303) 268-5407
 Toll free: +1 (866) 634-7464
 info@meisingerusa.com • www.meisingerusa.com

85FL139EN - 0321

