


# Twist Finishers

Finishing carbide burs with a twisted flute design, providing a remarkably smooth finish

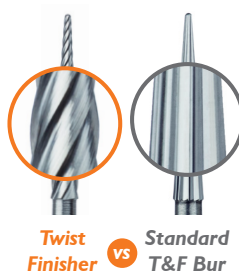

## Spiral Fluted Design

The spiral flute design allows the bur to stay in continuous contact with the tooth surface, which eliminates chatter caused by a standard T&F straight flute design.



## 12, 20, & 30 Flutes

Each size and shape is available in 12, 20, and 30 flutes to accommodate all esthetic trimming procedures.

MEISINGER's Twist Finisher Carbides are capable of producing a far smoother finish on composite restorations than is achievable with traditional straight-fluted finishing burs. The unique helical flute design ensures that the instrument stays in contact with the tooth surface when working on the restoration, eliminating chatter that occurs when using straight-fluted burs. This results in an extremely smooth finish with virtually no surface striations.

- Designed to complement all composites including hybrid composite materials
- Helical flute design ensures continuous contact, eliminating chatter
- Provides an extremely smooth surface with virtually no surface striations
- Made in Germany

- See a complete list of shapes and sizes on reverse side -

*Pair the Twist Finisher Carbides with the award winning Composite Twist Polishers for amazing composite restorations - See reverse side for details*

### MEISINGER USA

10150 E. Easter Avenue • Centennial, Colorado 80112 • USA  
Phone: +1 (303) 268-5400 • Fax: +1 (303) 268-5407  
customerservice@meisingerusa.com • www.meisingerusa.com

   @meisingerusa


*Meisinger* since 1888



# Composite Twist Polishers



For the best composite finish, pair the Twist Finisher Carbides with these Composite Twist Polishers.

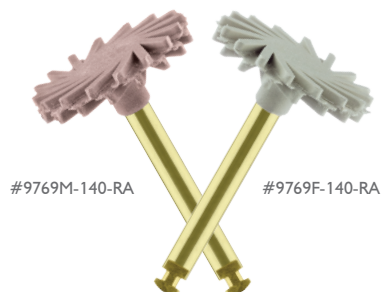

#9769M-140-RA

#9769F-140-RA

Contact your dealer representative to order



#LUS66

- Produces a brilliant natural luster
- The flexible shape design allows for easy and comfortable polishing on all surfaces
- Available in the LUS66 kit paired with the Twist Finisher Carbides (above)

NEW



## COMPOSITE SINGLES

Single-Use polishers designed to produce a brilliant luster on all composite materials.

- Quickly produces a brilliant luster on all composite materials
- Saves time from processing
- Unique shapes allow easy polishing on any surface

## Twist Finisher Sizes and Shapes

Fig.	HMD0132	HMD0132F	HMD0132U	HMD0133	HMD0133F	HMD0133U	HMD0134	HMD0134F	HMD0134U
Fine (12 Flute)	008			010			014		
Extra Fine (20 Flute)		008			010			014	
Ultra Fine (30 Flute)			008			010			014

Fig.	HMD0135	HMD0135F	HMD0135U	HMD0379	HMD0379F	HMD0379U	HMD0379	HMD0379F	HMD0379U
Fine (12 Flute)	014			018			023		
Extra Fine (20 Flute)		014			018			023	
Ultra Fine (30 Flute)			014			018			023

### MEISINGER USA

10150 E. Easter Avenue • Centennial, Colorado 80112 • USA  
 Phone: +1 (303) 268-5400 • Fax: +1 (303) 268-5407  
 customerservice@meisingerusa.com • www.meisingerusa.com

f t i @meisingerusa



Contact your local dealer representative to order today.