

Twist Polisher Kit for Composites

MEISINGER's Twist Polisher Kit has been specifically developed for easy high gloss polishing of all Composite restorations. The flexible polishing lamellae make it possible to adapt to any surface structure and also render flawless results in areas that are difficult to reach, such as occlusal surfaces. In addition, the two step polishers provide for short processing times and extremely high levels of durability.



2663

Twist Polishing Kit

Intraoral high-gloss polishing for all composite restorations

At a glance

- Flexible polishing lamellae can adapt to any surface structure
- Specifically developed for areas that are difficult to reach, such as occlusal surfaces
- Low heat production due to permanent air cooling of the material
- Designed for gentle action, which preserves the original surface structure
- Graduated diamond grit system allows for short processing times and a long service life
- User friendly and intuitive handling

MEISINGER USA, L.L.C.

10200 E. Easter Ave. • Centennial • Colorado • 80112

Tel.: +1 (303) 268-5400 • Fax: +1 (303) 268-5407

www.meisingerusa.com • customerservice@meisingerusa.com



twistfinishers *“ET” burs with a Twist*

MEISINGER's innovative Twist Finishers are capable of producing a far smoother finish on composite restorations than is achievable with traditional straight-fluted finishing burs. The unique helical design of the Twist Finishers ensures that the instrument stays in contact with the tooth surface when working on the restoration. This eliminates the chatter that occurs when using straight-fluted burs, resulting in an extremely smooth finish with virtually no surface striations left behind.

twistfinishers



Fig.	HMD0132	HMD0133	HMD0134	HMD0135	HMD0379
Diameter	008	010	014	014	023

In addition to all the shapes that are identical to the “ET” burs, MEISINGER also has a football shape!

Standard “ET” T&F Carbide Burs

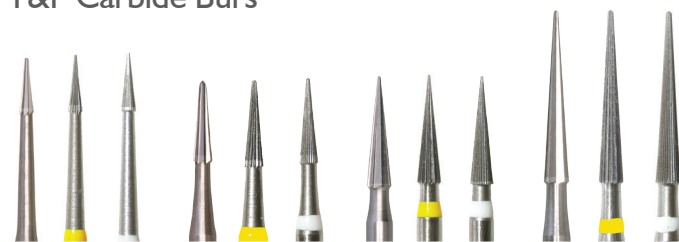


Fig.	HM132	HM133	HM134	HM135
Diameter	008	010	014	014


Available in 12, 20 and 30 flutes!

At a glance

- Helical design ensures better contact with surface during restoration process
- Extremely smooth finish with virtually no surface striations left behind
- Specifically designed for hybrid and new generation composite materials



VS.

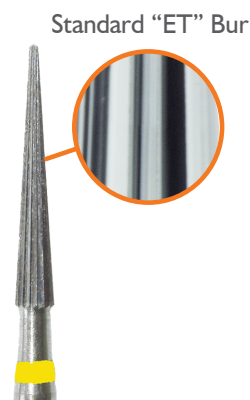

MEISINGER USA, L.L.C.

10200 E. Easter Ave. • Centennial • Colorado • 80112

Tel.: +1 (303) 268-5400 • Fax: +1 (303) 268-5407

www.meisingerusa.com • customerservice@meisingerusa.com

Meisinger
since
1888